









Introducing our most versatile and reliable cryo system

Building on the exceptional reliability of its predecessors, the Cryomatic MKII represents the perfect combination of reliability, automation, and clinical versatility. The innovative Cryomatic MK II features an intelligent coupling mechanism that seamlessly accommodates Keeler single-use and reusable probes without adapters, while the smart console automatically recognizes the probe type and self-adjusts, streamlining your workflow.

The intuitive interface features a backlit digital display and audio feedback system, providing clear, real-time information on freeze duration, gas cylinder status, and probe functionality. Designed for maximum clinical flexibility, the system is compatible with seven specialized reusable probes for various pathologies, plus a single-use probe that eliminates sterilization requirements. The Cryomatic MKII works efficiently with both CO2 and N2O medical gases, adapting to your specific needs and geographical requirements.

Focused on you: Local servicing and support

With over 50 year's experience in offering opthalmic cryosurgical units, Keeler has a long-standing legacy of trust and reliability when it comes to opthalmic cryosurgery.

We're proud to manufacture and assemble our instruments in the UK and the US, providing eyecare professionals with the highest quality ophthalmic devices, ensuring optical precision and reliability in every product. However, should you need assistance, our local expert technicians offer comprehensive care including repairs, annual servicing, and warranty support as well as training on the maintenance of your cryo systems.

Discover the benefits of the Cryomatic MK II



Compatibility with a range of Keeler probes

The console is compatible with Keeler's own range of reusable and single-use probes to provide users with the versatility to treat a range of pathologies.



Simple and convenient set up

The Cryomatic MKII is a complete solution for ophthalmic cryosurgery, featuring a simple graphic display, plus automated elements that allow users to be set up within a matter of minutes.



Automation for consistent performance

Compatible with Keeler's reusable and single-use probes, the integral software automatically detects and safely adjusts to your type of probe, applying the correct pressure levels to enable optimum freezing of the pathology.



Convenient, short purge cycle

The console ensures a short purge cycle, every time, to help to streamline your workflow, reducing preparation time between procedures. For added convenience, Keeler's single-use probes require no purge cycle at all, allowing for immediate use.



Purging facility with reusable probes

The purging facility removes any moisture or contaminants that could affect performance, ensuring consistent and reliable use throughout procedures.



Rapid freeze and quick defrost

The system delivers swift, controlled freezing, allowing users to precisely target pathological tissue with minimal impact to surrounding areas. This rapid response is matched by its quick defrost cycle, enabling smooth transitions between applications.

Probes that add versatility to your cryosurgery

The Cryomatic MK II is compatible with keeler's range of reusable and single-use probes. The Cryomatic MKII console and probes can be considered a superior treatment than using laser retinopexy for the following cases:

- · Treating peripheral retinal pathologies
- One-step procedure for retinal detachments. Laser retinopexy on the other hand, must be performed 24-48 hours after gas injection.
- Where there are obstructions within the eye limiting penetration of the laser, such as a cataract or vitreous opacity.



Keeler reusable probes

Keeler's range of reusable probes provides users with versatility to treat a range of pathologies. Designed for ease of connection, there are 8 tips to choose from depending on your requirements.

- A Standard retinal probe
- B Extended retinal probe
- C Mid-reach retinal probe
- Glaucoma probe
- E Curved cataract probe
- F Intra-vitreal probe
- G Collins trichiasis probe
- H Compact retinal probe

Keeler single-use probes

Keeler's single-use probes have been designed in collaboration with clinicians to provide users with optimal performance and safety when treating a range of pathologies such as retinal tears,

holes, Retinopathy of Prematurity (ROP), and pneumatic retinopexy. They can also be used during scleral buckle procedures for retinal detachments.

Keeler single-use probes offer the same quality, performance and precision as our renowned reusable retinal probes, whilst offering the peace of mind that a sterilised probe is on hand for every patient.



Orientation pip

Tactile feedback to ensure perfect placement of the tip.



Curved tip

Designed with the optimum angle to ensure accurate contact with the sclera.



Ergonomic design

The single-use probe is shaped to allow better control of the probe, with an insulated body for comfort during the freeze stage.





Increased safety

Constructed with a fail safe auto defrost feature for 100% confidence. No risk of forceful or sudden disassembly at any point in your cryosurgery procedure, ensuring that you and your patients are in safe hands.

Technical specifications


Reusable probes

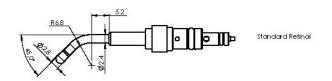

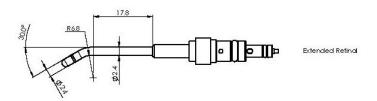


Standard retinal probe

SKU: 2509-P-8020

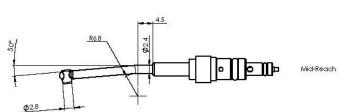

Extended retinal probe


SKU: 2509-P-8021

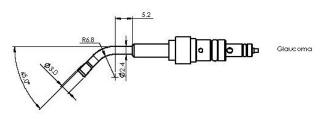

Mid-reach retinal probe

SKU: 2509-P-8022




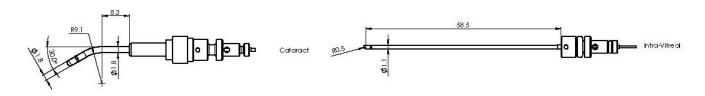


Curved cataract probe

SKU: 2509-P-8024



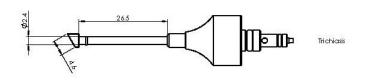

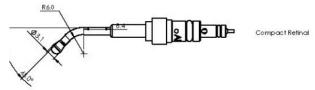

Collins trichiasis probe

SKU: 2509-P-8026

Compact retinal probe



Technical specifications

Cryomatic MK II

Cryo system	2509-P-1010
-------------	-------------

Gas specification	Medical grade nitrous oxide (N20) or medical grade carbon dioxide (CO2) in non-syphonic cylinders.
Operating range	3100-4480kPa (450-650PSI/31-45Bar)
Maximum cylinder pressure	8275kPa (1200PSI/83Bar)
Display	LCD Display (95mm x 70mm)

Electrical ratings

Input voltage range	100-240Vac (50/60Hz)
Power rating	15-30VA
Fuses	2 x T2AH 250V

Dimensions

Width	350mm (14")
Depth	200mm (8")
Height	190mm (7.5")
Weight	4.5kg (10lbs)

Classification and safety standards

Complies with	EN60601-1, UL60601-1, and CAN/CSA-C22.2 No 601-1
Equipment classification	Class 1, type BF (Applied Part)
Mode of operation	Continuous
Protection against ingress	Console IPx0 Footswitch IPX7

Technical specifications

Cylinder adapters

Your Cryomatic can be used with either Carbon Dioxide (CO²) or Nitrous Oxide (N²0) gas from non-syphonic bottles via one of the Pin Index Yoke (PIY) or screw thread adapters or cyclinder adapters.



SKU: 2508-P-7017

Pin Index Yoke for E Size Nitrous Oxide Gas Cylinder

SKU: 2508-P-7015

Pin Index Yoke for E Size Carbon Dioxide Gas Cylinder



SKU: 2508-P-7016

Adapter for F Size Carbon Dioxide Gas Cylinder

SKU: 2508-P-7018

Adapter for F Size Nitrous Oxide Gas Cylinder



Size E Cylinder

Dimensions (cm) 50x15 Gr. Weight (kg) 7

Size F Cylinder

Dimensions (cm) 86x14 Gr. Weight (kg) 18

Cryo accessories



SKU: 3304-K-1086

Cryo cart

What's in the box?

- · Cryomatic MKII console
- Footswitch
- Mains/power cord
- High-pressure gas hose
- Exhaust hose
- · Adjustable wrench
- · Instructions for use
- · 2x spare mains fuses



Visit keelerglobal.com for our full range including:













Contact us

Phone: 1-800-523-5620

Email: customerservice@keelerusa.com

UK CA 0120 C€ 1639

EC REP Keeler Europe Distribution, S.L. Colom, 453, nau D50 08223 Terrassa, Spain

СН REP



Keeler Ltd. Clewer Hill Road, Windsor, SL4 4AA



20022025 REV A